

Introduction

Prestained Protein Marker consists of high-purity and pre-stained recombinant proteins, which is convenient for monitoring of electrophoresis and transfer efficiency. This ready-to-load protein ladder is formulated in Laemmli buffer with 0.1% SDS, requiring no denaturation or reduction steps prior to SDS-PAGE or western blot applications.

Package Information

Components	C0115
Prestained Protein Marker	4×250 µl

Storage


Store at -20°C

Protocol

1. Thaw the marker to room temperature and mix thoroughly before use.
2. In a conventional mini or midi electrophoresis system, the loading volume of Prestained Protein Marker is 5-10 µl.
3. Store at -20°C immediately after use.




Note

1. This product is ready-to-use, no need to add reducing agent.
2. This product cannot be heated, which will cause the protein bands to degrade or discolor.
3. For your safety or health, please wear safety glasses, gloves, or protective clothing.

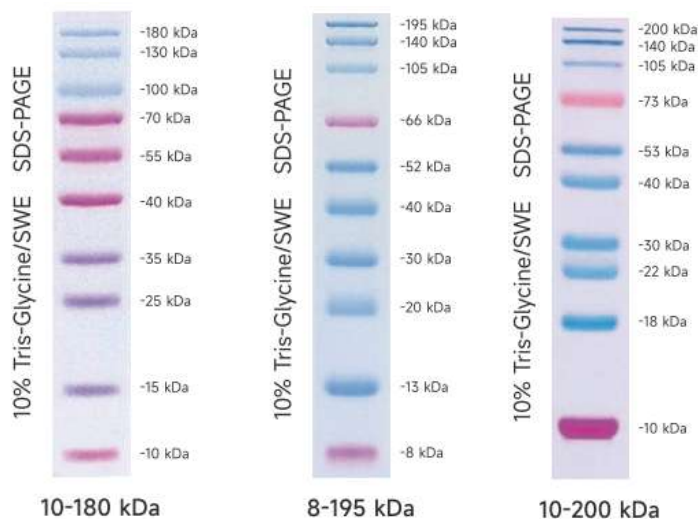

Prestained Protein Marker

Cat. #: C0115 Size: 4×250 µl

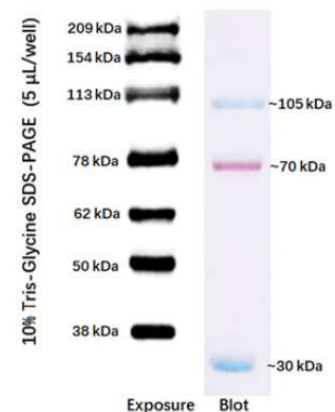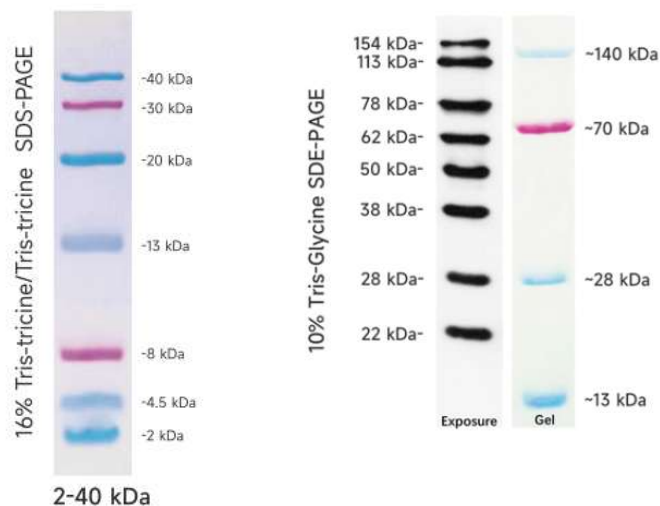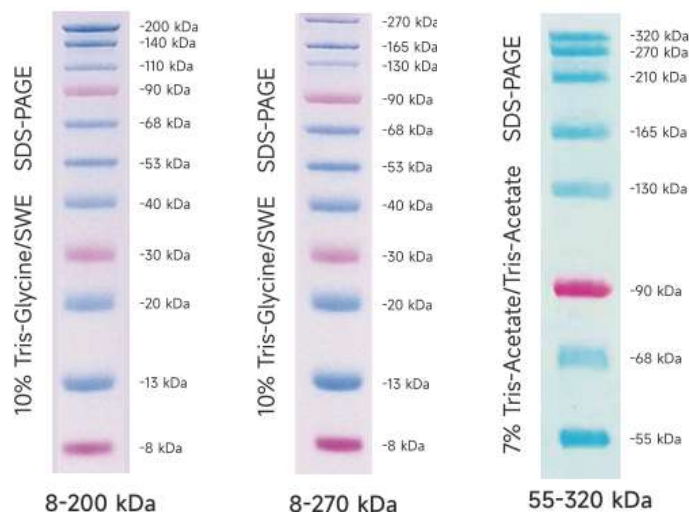

FOR RESEARCH USE ONLY